

Pro-Sport Devices Protocol: Knee Pain



Safety Precautions

Read all safety instructions in the owner's manual before use to prevent injury or death from improper use.

Intended Population: Is it for me?

- Read and understand English
- Have chronic (long-term) pain
- Have acute pain
- Have intractable (difficult to manage) pain
- Have post-surgical pain
- Have post-traumatic pain
- Not intended for those that have implanted electronic devices (demand-type cardiac pacemaker, implanted defibrillator, etc.)
- Not intended for those that have suspected or diagnosed heart conditions
- Not intended for those with suspected or diagnosed epilepsy
- Not intended for use when nursing, pregnant, or trying to become pregnant
- Not intended for use with undiagnosed pain syndromes
- Not intended for children (keep away from children)

Stop using the device if:

- Skin irritation occurs
- Hypersensitivity occurs
- Burning sensation occurs
- Feeling a headache
- Feeling lightheaded (people with low blood pressure should use with caution particularly in neck and lower abdomen)
- Feeling dizzy (people with low blood pressure should use with caution particularly in neck and lower abdomen)
- Fainting occurs (people with low blood pressure should use with caution particularly in neck and lower abdomen)
- Under the effects of drug or alcohol

Avazzia Device Contraindications

- Do not use on persons with an implanted pacemaker or other implanted electronic device.
- Do not use if pregnant or nursing.

Where can the device be used?

- Transcutaneously (through the skin) over normal, intact, clean, and healthy skin
- Sore and aching areas due to strain from exercise or normal household work activities
- Shoulder
- Waist
- Back
- Back of the neck
- Upper extremities (arms)
- Lower extremities (legs)
- Not intended for use over the thyroid or the carotid sinus nerves
- Not intended for use on the front of neck
- Not intended for use in the mouth or other internal orifices
- Not intended for use on the eyeballs
- Not intended for use transcranially (through the head)
- Not intended for use over menstruating or pregnant uterus
- Not intended for use over open wounds or rashes
- Not intended for use transthoracically (across or through the chest)
- Not intended for use over or in proximity to cancerous lesions
- Not intended for use over metallic implants present in the spine

Consult with your physician and use caution if:

- There is a tendency to hemorrhage following acute trauma or fracture
- Following recent surgical procedures when muscle contraction may disrupt the healing process
- Using over areas of the skin which lack normal sensation

Pro-Sport Devices


Protocol: Knee Pain



1. PLACE THE PADS

Attach the pads to the lead wire. Place one electrode pad on the inside of the knee at the joint line. Place another electrode pad on the outside of the knee at the joint line.

2 Pads

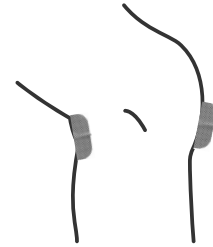

2. TURN ON THE DEVICE

Move the button (on the left side) to the "up" position.



3. SELECT THE "MODE"

Press the "up" or "down" arrow button to choose the mode. Press the "mode" button to confirm selection.


OR

4 Pads



4. INCREASE POWER LEVEL AND TREAT

Increase the power level by pressing the "+" key. Treat for 20 minutes.

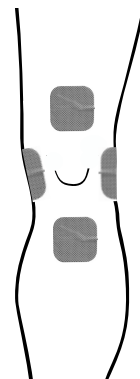

5. TURN OFF THE DEVICE REMOVE THE PADS RECORD TREATMENT NOTES



REMAIN SEATED OR LYING DOWN DURING TREATMENT. DO NOT STAND OR WALK.



KEEP DEVICE DRY